



Life Sciences

USD 2562b

Fluorodyne® II Filters

Bioburden, Sterility and Mycoplasma Control



Filtration. Separation. Solution.SM

Fluorodyne II Filter Series

Fluorodyne II filters feature a unique hydrophilic modified polyvinylidene fluoride (PVDF) membrane with low binding properties, broad chemical and temperature resistance and high flow rates.

Fluorodyne II Grade DBL Filters

These 0.45 micron rated filters are predominantly used in pre-filtration applications or viscous fluid filtrations. Validated removal of *Serratia marcescens* (ATCC[®] 14756) with a typical titer reduction (T_R) of $> 10^6$ ensures they are also suitable for bioburden control. A serial layer (0.65 / 0.45 micron) membrane construction allows for high throughputs and high flow rates compared with other 0.45 micron rated filters.




Fluorodyne II Grade DFL Filters

This validated 0.2 micron rated sterilizing grade filter retains *Brevundimonas diminuta* (ATCC[®] 19146) at 10^7 CFU/cm² Effective Filter Area (EFA) and permits very low extractables and adsorption. It is recommended for sterilizing filtration of pharmaceutical fluids including solutions with active ingredients, biologicals, biotech proteins, ophthalmics and other dilute preservative solutions.

Fluorodyne II Grade DJL Filters

This 0.1 micron rated filter with serial layer (0.2 / 0.1 micron) membrane construction assures high flow rates compared with other 0.1 micron (and even some 0.2 micron) rated filters. Fluorodyne II grade DJL filters are validated for removal of *Acholeplasma laidlawii* mycoplasma (ATCC[®] 23206) with a typical $T_R > 10^{11}$ and retention of *B. diminuta* (ATCC[®] 19146) at 10^7 CFU/cm² EFA. This allows for enhanced sterilization assurance as well as efficient mycoplasma control.

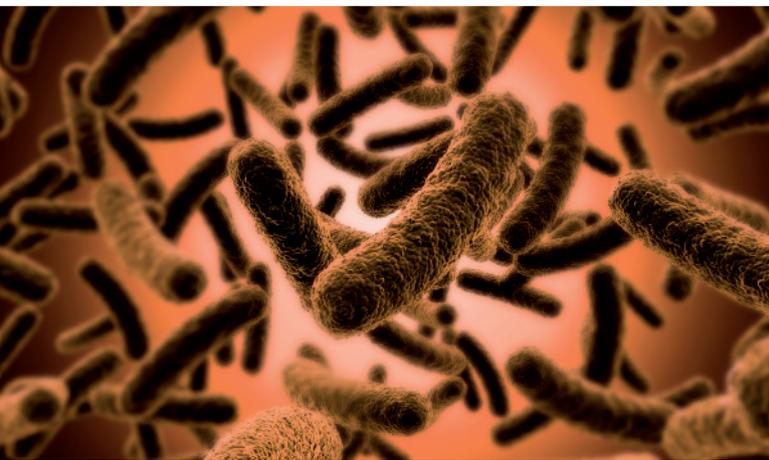





Allegro™ Systems: The Single-Use Solution



Cross-contamination elimination, assured sterility, reduction of manufacturing time, costs, and greater flexibility are clear objectives for the biopharmaceutical industry. These factors, coupled with increasing titers in drug manufacturing, may demand a new approach.

Allegro single-use systems featuring Kleepak filter capsules eliminate the need for cleaning and associated validation efforts, minimize major capital investments, increase flexibility and provide high assurance of product safety.

We provide full support for our single-use systems including training and validation services, to facilitate their use from upstream bioreactor, to final formulation and filling.

With a comprehensive range of scalable products in the Fluorodyne II filter series and Pall Allegro product platform, single-use systems incorporating Fluorodyne II filters can be used to process several milliliters up to large production-scale volumes.

Save Time, Get Results

Fluorodyne II filters are available in a wide range of scalable, encapsulated formats that allow for fast and easy scale-up, helping you rapidly deliver your products to the market.

Same Materials*

From laboratory-scale filters to production-scale assemblies, all Fluorodyne II filter products incorporate the same membrane and identical materials of construction, eliminating the need to requalify filter units as processes are scaled up.



Quality

Every Fluorodyne II pleated filter is:

- ▶ Integrity tested during manufacture
- ▶ Identified by lot and serial number for total traceability
- ▶ Supplied with a Certificate of Test confirming each filter:
 - Meets USP Biological Reactivity Test *in vivo*, for class VI-121 °C plastics
 - Meets cleanliness per USP Particulates in Injectables after flushing
 - Is non-fiber-releasing
 - Is non-pyrogenic per USP endotoxins (< 0.25 EU/mL)
 - Meets Total Organic Carbon (TOC) and water conductivity per USP purified water after flushing

* Except Novasip™ capsules. See material of construction tables for each product for further details.

	Mini Kleepak 20 Capsules Page 6
	Mini Kleepak Capsules Page 7
	Kleepak Capsules Page 8 - 9
	Novasip Capsules Page 10
	Junior Filter Cartridges MCY and SBF Style Page 11
	Junior Filter Cartridges Sealkleen Style Page 12
	High Area Filter Cartridges Page 13
	High Area Kleepak Nova Capsules Page 14 - 15



Mini Kleenpak 20 Capsules

Materials of Construction

Filter Membrane	Hydrophilic modified PVDF
Housing, Vent Plug and Support Material	Polypropylene
Filling Bell	Polycarbonate
Sealing Technology	Thermal bonding without adhesives

Operating Parameters⁽¹⁾

Maximum Operating Temperature and Pressure	1.4 bar (20 psi) at 20 °C
--	---------------------------

⁽¹⁾ In compatible fluids which do not soften, swell or adversely affect the filter or its materials of construction

Sterilization⁽²⁾

Pre-sterilized, subject to a minimum of 25 kGy of gamma irradiation

- ⁽²⁾
- Pre-sterilized Mini Kleenpak 20 filters must not be re-sterilized
 - Mini Kleenpak filters must not be sterilized in-situ by passing steam under pressure

Typical Hold up Volume

< 2.5 mL

Nominal Dimensions

Capsule Length	83 mm (3.3 in.)
Capsule Diameter	67 mm (2.7 in.)

Ordering Information⁽³⁾

Pall Part Number: KM5



P



S


Code	Filter Media
DFL	0.2 micron sterilizing grade
DJL	0.1 micron sterilizing grade and mycoplasma control

Inlet/Outlet Connection
 ¼ - ½ in. (6 - 13 mm) stepped hose barb with inner bore to accept female slip luer interior and outer diameter to accept filling bell outlet

Shipping Format
 Pre-sterilized using gamma irradiation

⁽³⁾ 3 filters per box

Typical Liquid Flow vs. Differential Pressure

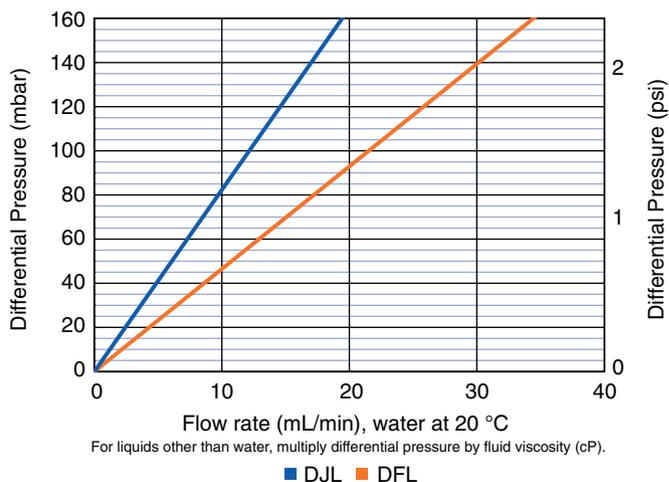

Nominal Effective Filter Area (EFA)

20 cm² (3.1 in.²)

Mini Kleenpak Capsules



Materials of Construction

Filter Membrane	Hydrophilic modified PVDF
Support/Drainage	Polypropylene
Capsule Shell	Polypropylene
Filling Bell	Polycarbonate
Sealing Technology	Thermal bonding without adhesives

Operating Parameters⁽¹⁾

Maximum Temperature	40 °C
Maximum Operating Pressure	4.1 bar (60 psi) at 40 °C
Maximum Differential Pressure (forward direction)	4.1 bar (60 psi) at 40 °C

⁽¹⁾ In compatible fluids which do not soften, swell or adversely affect the filter or its materials of construction

Sterilization⁽²⁾

Autoclave (G option only)	3 x 60 minutes at 125 °C
Gamma Irradiation (G option only)	Maximum of 50 kGy

- ⁽²⁾
- Pre-sterilized Mini Kleenpak capsules must not be re-sterilized.
 - Mini Kleenpak capsules must not be sterilized in-situ by passing steam under pressure

Typical Extractables in Water at 20 °C⁽³⁾

G option only	< 1 mg
S option only	< 5 mg

⁽³⁾ Tested on capsules without pre-flushing

Ordering Information⁽⁴⁾

Pall Part Number: KA02

Code	Filter Media
DFL	0.2 micron sterilizing grade
DJL	0.1 micron sterilizing grade and mycoplasma control


Code	Inlet/Outlet Connections
2	¼ - ½ in. (6 - 13 mm) stepped hose barb
8	½ - ¾ in. (13 - 19 mm) sanitary flange

Code	Shipping Format
G	Non-sterile gamma irradiatable/ autoclavable
S ⁽⁵⁾	Pre-sterilized using gamma irradiation

⁽⁴⁾ 3 filters per box

⁽⁵⁾ S option with Code 2 connection is provided with filling bell on outlet. It is removable for in-line use

Typical Liquid Flow vs. Differential Pressure

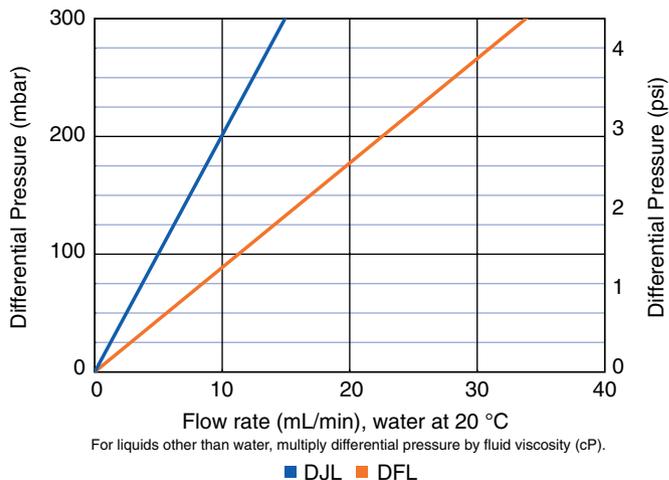

Nominal Dimensions

Maximum Diameter Including Valves	41 mm (1.6 in.)
Length - Code 2	105 mm (4.1 in.)
Length - Code 8	73 mm (2.9 in.)

Nominal Effective Filter Area (EFA)

200 cm² (0.22 ft²)

Kleenpak Capsules

Materials of Construction

Filter Membrane	Hydrophilic PVDF
Support/Drainage	Polypropylene
End Cap, Core and Cage	Polypropylene
Capsule Shell	Polypropylene
Sealing Technology	Thermal bonding without adhesives

Operating Parameters⁽¹⁾

Maximum Temperature	40 °C
Maximum Operating Pressure	5.2 bar (75 psi) at 20 °C 4.0 bar (58 psi) at 40 °C
Maximum Differential Pressure (forward direction)	4.0 bar (58 psi) at 40 °C

⁽¹⁾ In compatible fluids which do not soften, swell or adversely affect the filter or its materials of construction

Sterilization⁽²⁾

Autoclave (G and blank option only)	30 x 60 minutes at 135 °C 10 x 60 minutes at 140 °C
Gamma Irradiation (G option only)	Maximum of 50 kGy

- ⁽²⁾
- Pre-sterilized Kleenpak capsules must not be re-sterilized.
 - Kleenpak capsules must not be sterilized in-situ by passing steam under pressure

Typical Extractables in Water at 20 °C⁽³⁾

All Styles	< 5 mg per capsule
------------	--------------------

⁽³⁾ Tested on capsules without pre-flushing

Nominal Effective Filter Area (EFA)

KA1	400 cm ² (0.4 ft ²)
KA2	800 cm ² (0.9 ft ²)
KA3	1500 cm ² (1.6 ft ²)
KA4	3300 cm ² (3.6 ft ²)

Nominal Dimensions

	KA1	KA2	KA3	KA4
Diameter incl. Valves	94 mm (3.7 in.)	94 mm (3.7 in.)	105 mm (4.1 in.)	105 mm (4.1 in.)
Length - Code 1	117 mm (4.6 in.)	157 mm (6.2 in.)	174 mm (6.8 in.)	287 mm (11.3 in.)
Length - Code 2	158 mm (6.2 in.)	198 mm (7.8 in.)	—	—
Length - Code 6	157 mm (6.2 in.)	197 mm (7.7 in.)	210 mm (8.3 in.)	327 mm (12.8 in.)
Length - Code 12	137 mm (5.4 in.)	177 mm (7.0 in.)	—	—
Length - Code 16	137 mm (5.4 in.)	177 mm (7.0 in.)	192 mm (7.6 in.)	305 mm (12.0 in.)

Ordering Information

Pall Part Number:

Code	Code	Filter Media
1	DBL	0.45 micron bioburden control
2	DFL	0.2 micron sterilizing grade
3	DJL	0.1 micron sterilizing grade and mycoplasma control
4		


Code	Inlet/Outlet Connections
1	1/2 in. (25 - 38 mm) sanitary flange
2	1/4 - 1/2 in. (6 - 13 mm) stepped hose barb
6	1/2 in. (13 mm) single hose barb
16	1/2 in. sanitary flange inlet and 1/2 in. (13 mm) single hose barb outlet

Code	Shipping Format
G	Non-sterile gamma irradiatable/autoclavable
S	Pre-sterilized using gamma irradiation
Blank ⁽⁴⁾	Non-sterile autoclavable


⁽⁴⁾ Only available with DBL and DFL filter media




Typical Liquid Flow vs. Differential Pressure (KA1)

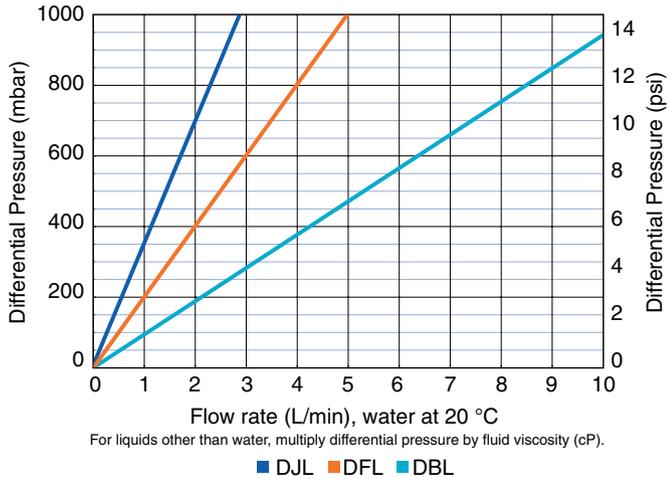


Typical Liquid Flow vs. Differential Pressure (KA2)

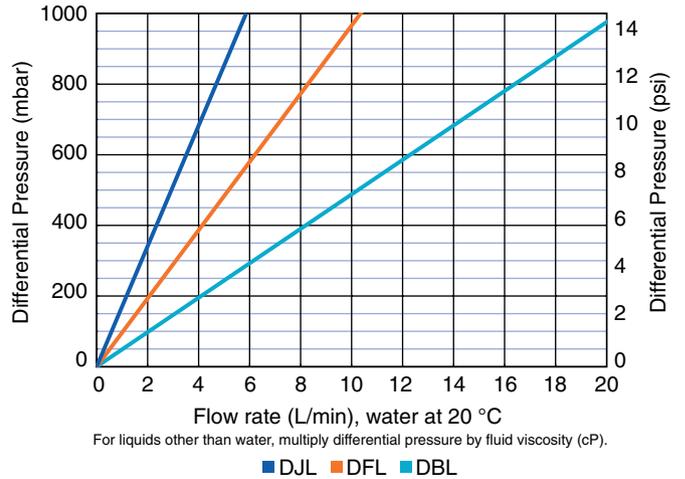

Typical Liquid Flow vs. Differential Pressure (KA3)

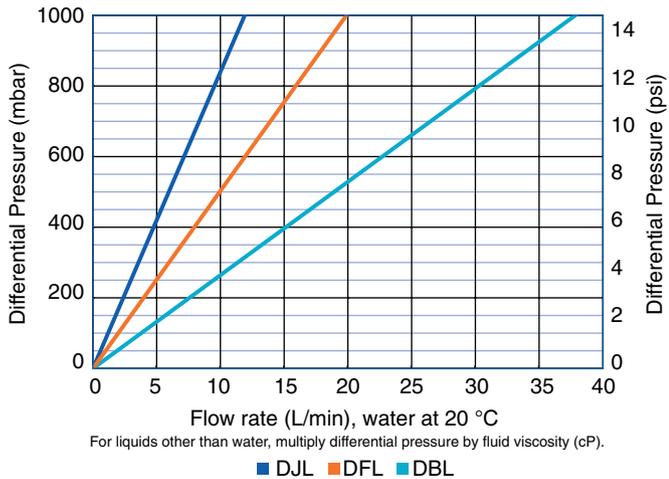

Typical Liquid Flow vs. Differential Pressure (KA4)

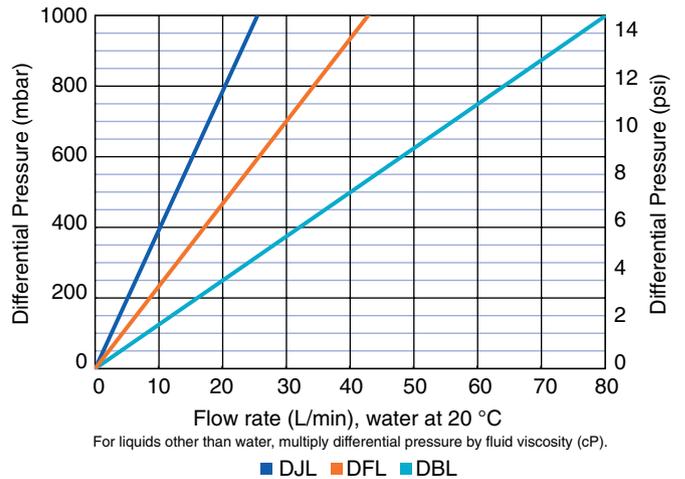

Novasip Capsules



Materials of Construction

Filter Membrane	Hydrophilic modified PVDF
Support/Drainage	Polypropylene
End Cap, Core and Cage	Polypropylene
Capsule Bowl	Polyetherimide
Sealing Technology	Thermal bonding without adhesives
Housing Head	Polyetherimide with TiO ₂ whitener ⁽¹⁾

⁽¹⁾ TiO₂ is an insoluble inorganic mineral filler that does not contribute to organic extractables

Operating Parameters⁽²⁾

Maximum Temperature	60 °C
Maximum Operating Pressure	6.5 bar (94 psi) at 40 °C 2.0 bar (29 psi) at 60 °C
Maximum Differential Pressure (forward direction)	5.3 bar (77 psi) at 40 °C 4.1 bar (60 psi) at 60 °C

⁽²⁾ In compatible fluids which do not soften, swell or adversely affect the filter or its materials of construction

Sterilization

Autoclave or <i>In Situ</i> Steam	30 x 60 minutes at 125 °C 10 x 60 minutes at 140 °C
-----------------------------------	--

Typical Extractables in Water at 20 °C⁽³⁾

< 5 mg per capsule

⁽³⁾ Tested on capsules without pre-flushing

Ordering Information

Pall Part Number:


C

P

⁽⁴⁾ DJL filter media is not available in Size 1

Code	Holder Volume	Size	Code	Filter Media	Connection Options	Code	Vent/Drain
L	Special insert for low hold up volume	05	DFL	0.2 micron sterilizing grade	1 ½ in. (25 - 38 mm) sanitary flange	Blank	Vent: Quick connect and disconnect coupling (compatible with Stäubli* fitting) Drain: Hose barb for ½ in. - ¼ in. (4 - 6 mm) i.D. tube, with valve
LM	Small Bowl (Only for size 05)	1	DJL ⁽⁴⁾	0.1 micron sterilizing grade and mycoplasma control		A	Vent and drain: Quick connect and disconnect coupling (Stäubli compatible) with valve
Blank	Standard Capsule	3				B	Vent and drain: ½ in. (13 mm) sanitary flange, no valve sanitary clamp

Typical Liquid Flow vs. Differential Pressure

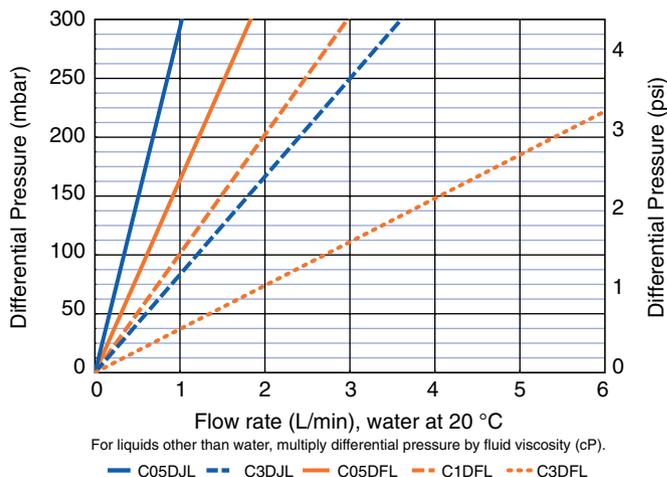

Nominal Dimensions

	CLM05	C05/C1/C3
Diameter Incl. Valves	123 mm (4.8 in.)	123 mm (4.8 in.)
Overall Length	84 mm (3.3 in.)	157 mm (6.2 in.)

Nominal Effective Filter Area (EFA)

C(LM)05	400 cm ² (0.4 ft ²)
C(L)1	700 cm ² (0.7 ft ²)
C(L)3	1500 cm ² (1.6 ft ²)

Junior Filter Cartridges - MCY and SBF Style



Materials of Construction

Filter Membrane	Hydrophilic modified PVDF
Support/Drainage	Polypropylene
End Cap, Core and Cage	Polypropylene
Sealing Technology	Thermal bonding without adhesives
O-rings	Silicone elastomer

Operating Parameters⁽¹⁾

Maximum Temperature	80 °C
Maximum Differential Pressure (forward direction)	5.3 bar (77 psi) at 40 °C 3.4 bar (50 psi) at 80 °C

- ⁽¹⁾ • In compatible fluids which do not soften, swell or adversely affect the filter or its materials of construction
• DFL filter media may be used at 90 °C. Contact Pall for further details

Sterilization

Autoclave and <i>In Situ</i> Steam	30 x 60 minutes at 125 °C 10 x 60 minutes at 140 °C ⁽²⁾
------------------------------------	---

⁽²⁾ MCY style only

Typical Extractables in Water at 20 °C⁽³⁾

< 5 mg per Junior cartridge

⁽³⁾ Tested on elements without pre-flushing

Nominal Dimensions


	MCY1110	MCY2220	MCY4440	SBF
Length	32 mm (1.3 in.)	57 mm (2.2 in.)	106 mm (4.2 in.)	51 mm (2.0 in.)

Ordering Information

*Pall Part Number:

Code	Style
MCY1110	Single open-ended with double external o-rings and locating lugs
MCY2220	Single open-ended with double external o-rings and locating lugs
MCY4440	Single open-ended with double external o-rings and locating lugs
SBF ⁽⁴⁾	Single open-ended with double external o-rings

Typical Liquid Flow vs. Differential Pressure

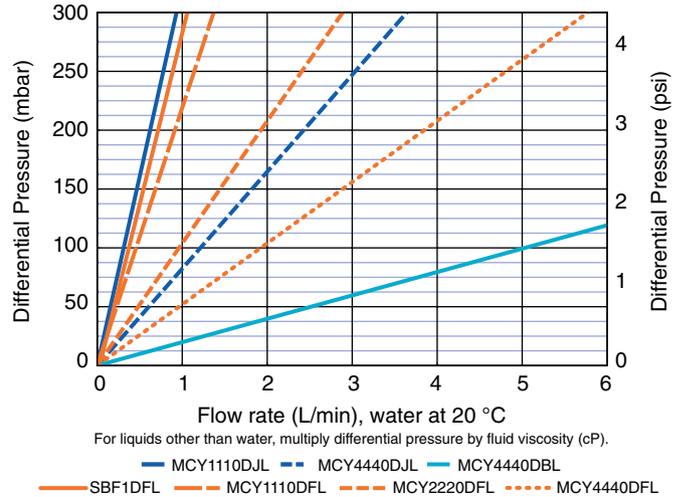

Nominal Effective Filter Area (EFA)

MCY1110	400 cm ² (0.4 ft ²)
MCY2220	700 cm ² (0.7 ft ²)
MCY4440	1500 cm ² (1.6 ft ²)
SBF	400 cm ² (0.4 ft ²)

P	H4
	O-ring Material
	Silicone elastomer (standard)
	⁽⁴⁾ Only available with DFL filter media
	⁽⁵⁾ Only available in MCY4440 style
	⁽⁶⁾ Only available in MCY1110 and MCY4440 styles

Code	Filter Media
DBL ⁽⁶⁾	0.45 micron bioburden control
DFL	0.2 micron sterilizing grade
DJL ⁽⁶⁾	0.1 micron sterilizing grade and mycoplasma control

Junior Filter Cartridges - Sealkleen Style



Materials of Construction

Filter Membrane	Hydrophilic modified PVDF
Support/Drainage	Polypropylene
End Cap, Core and Cage	Polypropylene
Sealing Technology	Thermal bonding without adhesives

Operating Parameters⁽¹⁾

Maximum Temperature	80 °C
Maximum Differential Pressure (forward direction)	5.3 bar (77 psi) at 50 °C 3.4 bar (50 psi) at 80 °C

- ⁽¹⁾
- In compatible fluids which do not soften, swell or adversely affect the filter or its materials of construction
 - DFL filter media may be used at 90 °C. Contact Pall for further details

Sterilization

Autoclave or <i>In Situ</i> Steam	30 x 60 minutes at 125 °C
-----------------------------------	---------------------------

Typical Extractables in Water at 20 °C⁽²⁾


< 5 mg (per Sealkleen cartridge)

⁽²⁾ Tested on elements without pre-flushing

Nominal Dimensions

	SLK7001	SLK7002
Overall Length	66 mm (2.6 in.)	133 mm (5.2 in.)

Typical Liquid Flow vs. Differential Pressure

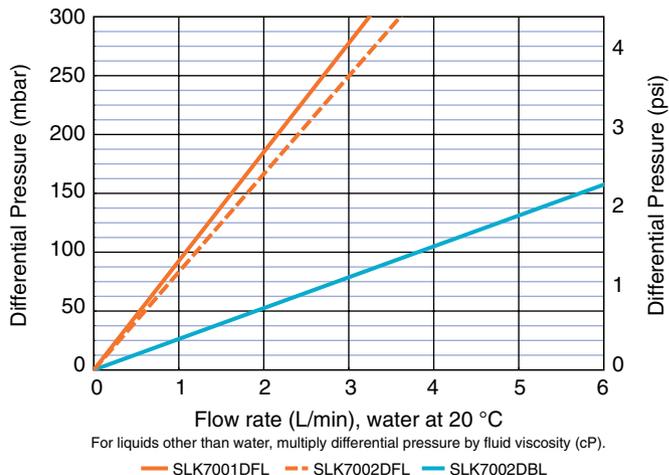

Nominal Effective Filter Area (EFA)

SLK7001	900 cm ² (1.0 ft ²)
SLK7002	1900 cm ² (2.0 ft ²)

Ordering Information

Pall Part Number: SLK P

Code	Code	Filter Media
7001 ⁽³⁾	DBL	0.45 micron bioburden control
7002	DFL	0.2 micron sterilizing grade

⁽³⁾ Only available with DFL filter media

High Area Filter Cartridges



Materials of Construction

Filter Membrane	Hydrophilic PVDF
Support/Drainage	Polypropylene
Core/End Caps	Polypropylene
Cage (DFL and DJL filter media)	Polypropylene with TiO ₂ whitener ⁽¹⁾
Cage (DBL filter media)	Polypropylene
O-rings	Silicone elastomer
Sealing Technology	Thermal bonding without adhesives

⁽¹⁾ TiO₂ is an insoluble inorganic mineral filler that does not contribute to organic extractables

Operating Parameters⁽²⁾

Maximum Temperature	80 °C
Maximum Differential Pressure (forward direction)	5.3 bar (77 psi) at 50 °C 3.4 bar (49 psi) at 80 °C
Maximum Differential Pressure (reverse direction)	0.3 bar (4 psi)

⁽²⁾ • In compatible fluids which do not soften, swell or adversely affect the filter or its materials of construction
• DFL filter media may be used at 90 °C. Contact Pall for further details

Sterilization

Autoclave or <i>In Situ</i> Steam	30 x 60 minutes at 125°C 10 x 60 minutes at 140°C
-----------------------------------	--

Typical Extractables in Water at 20 °C⁽³⁾

< 5 mg per 254 mm (10 in.) module

⁽³⁾ Tested on elements without pre-flushing

Ordering Information

Pall Part Number:

AB

Code	Nominal Length
05	125 mm (5 in.)
1	254 mm (10 in.)
2	508 mm (20 in.)
3	762 mm (30 in.)

Code	Filter Media
DBL	0.45 micron bioburden control
DFL	0.2 micron sterilizing grade
DJL	0.1 micron sterilizing grade and mycoplasma control


Code	Adapter Style
7	Pall Code 7 double O-ring bayonet lock and fin
2	Pall Code 2 double O-ring bayonet lock, no fin (Code 05 only)

P

H4

O-ring material
Silicone elastomer (Other materials available on request)

Typical Liquid Flow vs. Differential Pressure

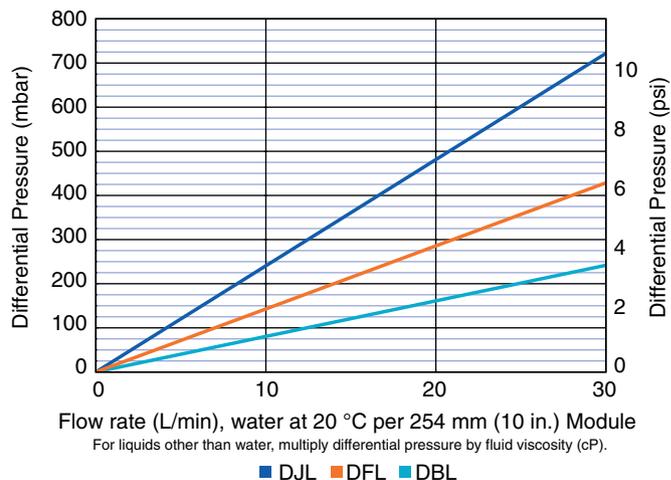

Integrity Test Values (Air test gas, water wet)

Max. allowable forward flow values for 254 mm (10 in.) module at 20°C	
DBL	13 mL/min at 1240 mbar (18 psi)
DFL	12 mL/min at 2760 mbar (40 psi)
DJL	29 mL/min at 4475 mbar (65 psi)

Contact Pall for multi-element integrity test values and recommended test procedures

Nominal Effective Filter Area (EFA) per 254 mm (10 in.) Module

DBL	5500 cm ² (5.9 ft ²)
DFL	5500 cm ² (5.9 ft ²)
DJL	5500 cm ² (5.9 ft ²)

High Area Kleenpak Nova Capsules

Materials of Construction

Filter Membrane	Hydrophilic modified PVDF
Support/Drainage	Polypropylene
Core/End Caps	Polypropylene
Cage (DFL and DJL filter media)	Polypropylene with TiO ₂ whitener ⁽¹⁾
Cage (DBL filter media)	Polypropylene
O-rings	Silicone elastomer
Sealing Technology	Thermal bonding without adhesives
Housing Bowl	Polypropylene
Housing Head	Polypropylene

⁽¹⁾ TiO₂ is an insoluble inorganic mineral filler that does not contribute to organic extractables

Operating Parameters⁽²⁾

Maximum Temperature	40 °C
Maximum Operating Pressure	3 bar (44 psi) at 40 °C
Maximum Differential Pressure (forward direction)	3 bar (44 psi) at 40 °C

⁽²⁾ In compatible fluids which do not soften, swell or adversely affect the filter or its materials of construction

Sterilization⁽³⁾

Autoclave (G option only)	3 x 60 minutes at 125 °C
Gamma Irradiation (G option only)	Maximum of 50 kGy

⁽³⁾ • Pre-sterilized Kleenpak Nova capsules must not be re-sterilized
• Kleenpak Nova capsules must not be sterilized in-situ by passing steam under pressure

Typical Extractables in Water at 20 °C⁽⁴⁾

G option	< 10 mg per 254 mm (10 in.) module
S option	< 20 mg per 254 mm (10 in.) module

⁽⁴⁾ Tested on capsules without pre-flushing

Integrity Test Values (Air test gas, water wet)⁽⁵⁾

Max. allowable forward flow values for 254 mm (10 in.) filter at 20°C	
DBL	13 mL/min at 1240 mbar (18 psi)
DFL	12 mL/min at 2760 mbar (40 psi)
DJL	29 mL/min at 4475 mbar (65 psi)

⁽⁵⁾ Contact Pall for multi-element integrity test values and recommended test procedures

Nominal Effective Filter Area (EFA) per 254 mm (10 in.) Module


DBL	5500 cm ² (5.9 ft ²)
DFL	5500 cm ² (5.9 ft ²)
DJL	5500 cm ² (5.9 ft ²)

Nominal Dimensions


In-line	NP5	NP6	NP7	NP8
Maximum diameter including valves	154 mm (6.1 in.)			
Length with hose barb inlet/outlet	275 mm (10.8 in.)	397 mm (15.6 in.)	644 mm (25.4 in.)	895 mm (35.2 in.)
Length with sanitary inlet/outlet	213 mm (8.4 in.)	335 mm (13.2 in.)	584 mm (23.0 in.)	834 mm (32.8 in.)
T-style	NT6	NT7	NT8	
Maximum diameter including valves	—	240 mm (9.5 in.)	240 mm (9.5 in.)	240 mm (9.5 in.)
Length	—	349 mm (13.7 in.)	598 mm (23.5 in.)	848 mm (33.4 in.)



Typical Liquid Flow vs. Differential Pressure (NP)

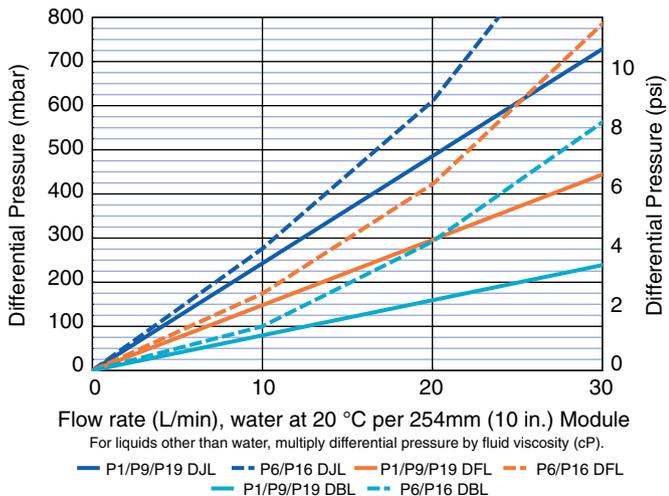

Typical Liquid Flow vs. Differential Pressure (NT)

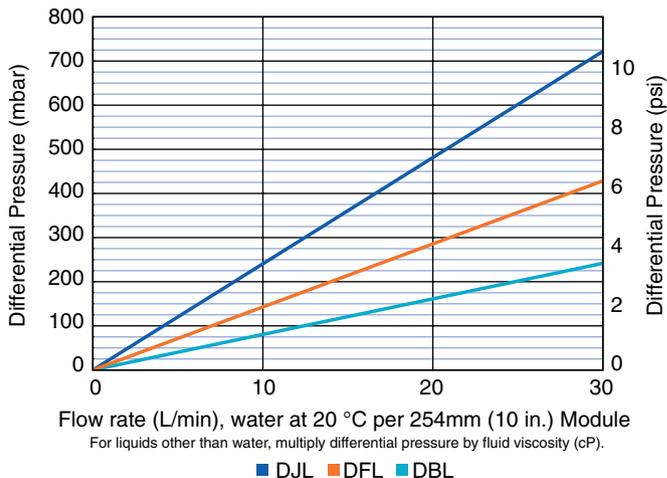

Ordering Information

Pall Part Number:

N

Code	Style
P	In-line
T	T-style

Code	Filter Size
5 ⁽¹⁾	127 mm (5 in.)
6	254 mm (10 in.)
7	508 mm (20 in.)
8	762 mm (30 in.)

Code	Filter Media
DBL	0.45 micron bioburden control
DFL	0.2 micron sterilizing grade
DJL	0.1 micron sterilizing grade and mycoplasma control

P

Code	Shipping Format
G	Non-sterile gamma irradiatable/ autoclavable
S	Pre-sterilized using gamma irradiation

⁽¹⁾ Only available in In-Line style and with DJL and DFL filter media only

⁽²⁾ For In-line (Code P) only

⁽³⁾ For T-style (Code T) only

Code	Connection Options
1	1 – 1 ½ in. (25 - 38 mm) sanitary flange inlet and outlet
9	1 in. (25 mm) single barb hose barb inlet and outlet
19	1 – 1 ½ in. (25 - 38 mm) sanitary flange inlet and 1 in. (25 mm) single barb hose barb outlet
6 ⁽²⁾	½ in. (13 mm) single barb hose barb inlet and outlet
16 ⁽²⁾	1 – 1 ½ in. (25 - 38 mm) sanitary flange inlet and ½ in. (13 mm) single barb hose barb outlet
1H ⁽³⁾	1 – 1 ½ in. (25 - 38 mm) sanitary flange inlet and outlet, with ½ in. sanitary port on inlet
1H9 ⁽³⁾	1 – 1 ½ in. sanitary flange inlet and 1 in. (25 mm) single barb hose barb outlet with ½ in. sanitary port on inlet

Code	Vent/Drain
Blank	Stäubli* vent and stepped hose barb drain
A	Stäubli vent and drain



Life Sciences

United States

800.717.7255 toll free (USA)
516.484.5400 phone
516.801.9548 fax
biopharm@pall.com E-mail

Europe

+41 (0)26 350 53 00 phone
+41 (0)26 350 53 53 fax
LifeSciences.EU@pall.com E-mail

Filtration. Separation. Solution.SM



Visit us on the Web at www.pall.com/biopharm

E-mail us at biopharm@pall.com

International Offices

Pall Corporation has offices and plants throughout the world in locations such as: Argentina, Australia, Austria, Belgium, Brazil, Canada, China, France, Germany, India, Indonesia, Ireland, Italy, Japan, Korea, Malaysia, Mexico, the Netherlands, New Zealand, Norway, Poland, Puerto Rico, Russia, Singapore, South Africa, Spain, Sweden, Switzerland, Taiwan, Thailand, the United Kingdom, the United States, and Venezuela. Distributors in all major industrial areas of the world.

The information provided in this literature was reviewed for accuracy at the time of publication. Product data may be subject to change without notice. For current information consult your local Pall distributor or contact Pall directly.

© 2009, Pall Corporation. Pall, , Allegro, the Allegro Design, Fluorodyne, Kleenpak and Novasip are trademarks of Pall Corporation. ® indicates a trademark registered in the USA and TM indicates a common law trademark.

Filtration. Separation. Solution.SM and UpScaleSM are service marks of Pall Corporation. *ATCC is a registered trademark of American Type Culture Collection. Stäubli is a trademark of Stäubli AG.

11/2017, UK GN17.06971

USD 2562b