

ÄKTA pure

Unpacking Instruction

Read this before unpacking ÄKTA pure



This instruction contains information that is important for safety during the unpacking of ÄKTA™ pure. Read the whole instruction carefully before starting to unpack the system.

When you receive the delivery

- Record on the receiving documents if there is any apparent damage on the delivery box. Inform your Cytiva representative of such damage.
- Move the delivery box to a protected location indoors.

About this instruction

This instruction covers unpacking of ÄKTA pure. See ÄKTA pure Operating instructions for instructions on how to install the instrument.

Instruction

Follow the instruction below to unpack the instrument and lift the instrument onto the bench.



CAUTION

Heavy object. Use proper lifting equipment, or use two or more persons when moving the instrument. All lifting and moving must be performed in accordance with local regulations.

Step	Action
------	--------

- | | |
|---|------------------------------------|
| 1 | Cut and remove the plastic straps. |
|---|------------------------------------|



Step	Action
------	--------

- | | |
|---|--|
| 2 | Lift off and remove the lid and protective foam. |
|---|--|



- | | |
|---|---|
| 3 | Check the contents in the Buffer tray, and lift off the packages from the tray. |
|---|---|

- | | |
|---|---|
| 4 | Lift off the cardboard hood and remove the protecting material from the instrument. |
|---|---|




Step	Action
------	--------

- | | |
|---|--|
| 5 | Lift off the tray on the wet side on the front of the instrument to access the instrument handles. |
|---|--|



Step	Action
------	--------

- | | |
|---|--|
| 6 | Prepare for lifting. Use two or more persons and grip the instrument from the front, from the back or from either side (only one side is shown below): |
|---|--|

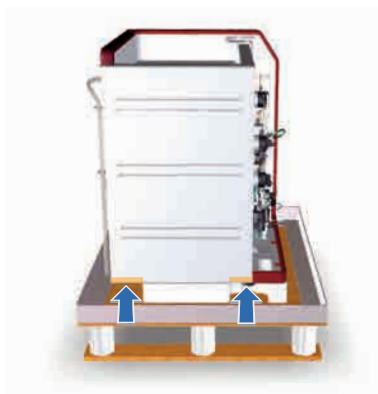

Step	Action
------	--------

- | | |
|---|---|
| 7 | Lift the instrument over the foam attached to the plywood board, and pull away the board from under the instrument. |
|---|---|



- | | |
|---|---|
| 8 | Dispose of the packaging material in accordance with local regulations. |
|---|---|

Page intentionally left blank



cytiva.com/akta

Cytiva and the Drop logo are trademarks of Global Life Sciences IP Holdco LLC or an affiliate.

AKTA is a trademark of Global Life Sciences Solutions USA LLC or an affiliate doing business as Cytiva.

All other third-party trademarks are the property of their respective owners.

© 2020 Cytiva

All goods and services are sold subject to the terms and conditions of sale of the supplying company operating within the Cytiva business. A copy of those terms and conditions is available on request. Contact your local Cytiva representative for the most current information.

For local office contact information, visit cytiva.com/contact

29020657 AD V:5 11/2020