

# Purification of antibodies using ÄKTA start and HiTrap Protein G HP column

**Intellectual Property Notice:** The Biopharma business of GE Healthcare was acquired by Danaher on 31 March 2020 and now operates under the Cytiva™ brand. Certain collateral materials (such as application notes, scientific posters, and white papers) were created prior to the Danaher acquisition and contain various GE owned trademarks and font designs. In order to maintain the familiarity of those materials for long-serving customers and to preserve the integrity of those scientific documents, those GE owned trademarks and font designs remain in place, it being specifically acknowledged by Danaher and the Cytiva business that GE owns such GE trademarks and font designs.

## **cytiva.com**

GE and the GE Monogram are trademarks of General Electric Company.  
Other trademarks listed as being owned by General Electric Company contained in materials that pre-date the Danaher acquisition and relate to products within Cytiva's portfolio are now trademarks of Global Life Sciences Solutions USA LLC or an affiliate doing business as Cytiva.  
Cytiva and the Drop logo are trademarks of Global Life Sciences IP Holdco LLC or an affiliate.  
All other third-party trademarks are the property of their respective owners.  
© 2020 Cytiva  
All goods and services are sold subject to the terms and conditions of sale of the supplying company operating within the Cytiva business. A copy of those terms and conditions is available on request. Contact your local Cytiva representative for the most current information.  
For local office contact information, visit [cytiva.com/contact](https://www.cytiva.com/contact)

# Purification of antibodies using ÄKTA™ start and HiTrap™ Protein G HP column

**Polyclonal and monoclonal antibodies (MAbs) from many species can be purified using HiTrap Protein G Sepharose™ affinity columns. In this Application note, Keisha uses ÄKTA start system and a HiTrap Protein G HP column to reduce the time taken for one-step purification of MAbs from culture supernatant to high purity (> 95%).**


Keisha switched to ÄKTA start system (comprising the ÄKTA start instrument, UNICORN™ start 1.0 control software and the Frac30 fraction collector) and selected a HiTrap Protein G HP column for purifying the MAb. HiTrap Protein G HP is a convenient, ready-to-use column prepacked with Protein G Sepharose High Performance chromatography medium (resin).

## Introduction

The diversity of antibody-antigen interactions has created many uses for antibodies. They are used for therapeutic and diagnostic applications as well as for immunochemical techniques within general research.

Keisha is a researcher at a start-up biotech company focusing on characterizing MAbs at a preclinical phase of development. Obtaining highly purified antibodies is of prime importance for her work. Keisha was often challenged by the manual purification methods for purifying MAbs. Too much effort was being spent on time-consuming work such as loading large sample volumes manually. As a result, she was able to purify only a single batch of MAb per day, and was obtaining inconsistent yields and varying degrees of purity across batches.

Keisha purifies MAbs in a single step using Protein G affinity chromatography. The basis for antibody affinity purification using protein G is the high affinity and specificity for the Fc region of IgG from a variety of species. Protein G media bind IgG from a broader range of eukaryotic species and bind more subclasses of IgG than Protein A. Protein G has been immobilized to several different matrixes, resulting in an excellent means of isolating IgG and IgG subclasses from ascites, cell culture supernatants, and serum.

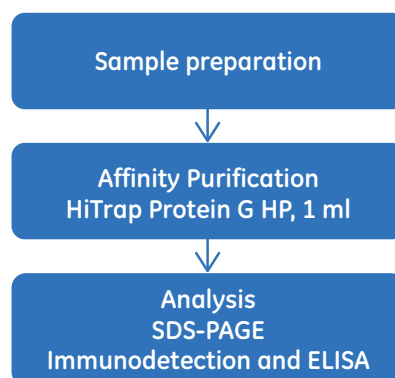

**Fig 1.** Workflow illustrating the steps involved in purifying MAbs.



Methods

Keisha followed the workflow described in Figure 1. The sample was prepared by diluting the culture supernatant two-fold using binding buffer (20 mM sodium phosphate, 150 mM NaCl, pH 7.0). The diluted sample was kept at room temperature for 15 to 20 min to ensure solubility of the protein.

Large sample volumes can be loaded on ÄKTA start using Superloop™.

Keisha used UNICORN start to modify the predefined affinity template, as described in Table 1, and used the System Control module of the software to run the method.

The ÄKTA start pump was used to load 20 ml of diluted culture supernatant onto a HiTrap Protein G HP 1 ml column. After the sample was loaded, the column was washed with 10 column volumes (CV) of binding buffer.

Create custom methods using drag-and-drop functionality in UNICORN start Method Editor.

The bound proteins were eluted by step gradient using elution buffer (100 mM glycine-HCl, pH 2.7). The first step was 5 CV of 70% B to remove nonspecific contaminants and the second step was 10 CV of 100% B to elute the protein of interest. Peak fractions were collected using Frac30 fraction collector and pooled following the purification run. The purity and specificity of the purified MAb was assessed by SDS-polyacrylamide gel electrophoresis (SDS-PAGE, using a 12.5% polyacrylamide gel) and immunodetection.

Table 1. UNICORN start method overview

Method flow	Method settings
Method settings	CV <sup>1</sup> Pressure limit 0.3 MPa <sup>2</sup> Flow rate 1 ml/min
Prime and equilibration	5 CV <sup>1</sup>
Sample application	Apply sample using pump Sample volume 20 ml
Wash out unbound	10 CV <sup>1</sup>
Elution and fractionation	Gradient elution: Step 1: 70% B 5 CV <sup>1</sup> Step 2: 100% B 10 CV <sup>1</sup>  Fractionation: Peak fractionation level based (start/end: 5 mAU)
Prime and equilibration	5 CV <sup>1</sup>

<sup>1</sup> CV = Column volumes  
<sup>2</sup> 0.3 MPa = 3 bar (43.5 psi)

Results

Using ÄKTA start with UNICORN start software and Frac30 fraction collector, Keisha was able to purify MAb from culture supernatant in a single purification step, and she could identify and pool the peak fractions (Fig 2A). ÄKTA start system with the affinity method gave high yields and > 95% purity of the target protein. SDS-PAGE and immunoblot analysis of the purified sample with protein-specific antibodies confirmed the purity and specificity of the MAb, with no contaminating protein bands (Fig 2B). ELISA experiments confirmed that the purified MAb reacted specifically with the antigen, and activity was detected at concentrations as low as 6.25 ng/ml (data not shown).

Ready-to-use methods are available for column preparation and system cleaning.

ÄKTA start system reduced the manual intervention required for MAb purification. The time saved allowed Keisha to conduct several purification runs per day, and she was able to achieve consistent results across runs.

Column: HiTrap Protein G HP, 1 ml  
Sample: Culture supernatant diluted in 2x binding buffer, pH 7  
Sample load: 20 ml  
Binding buffer: 20 mM sodium phosphate, 150 mM NaCl, pH 7.0  
Elution buffer: 100 mM glycine-HCl, pH 2.7  
Flow rate: 1 ml/min  
Gradient: Step 1: 70% B, 5CV  
Step 2: 100% B, 10 CV  
System: ÄKTA start, UNICORN start and Frac30 fraction collector  
Detection: UV (280 nm)

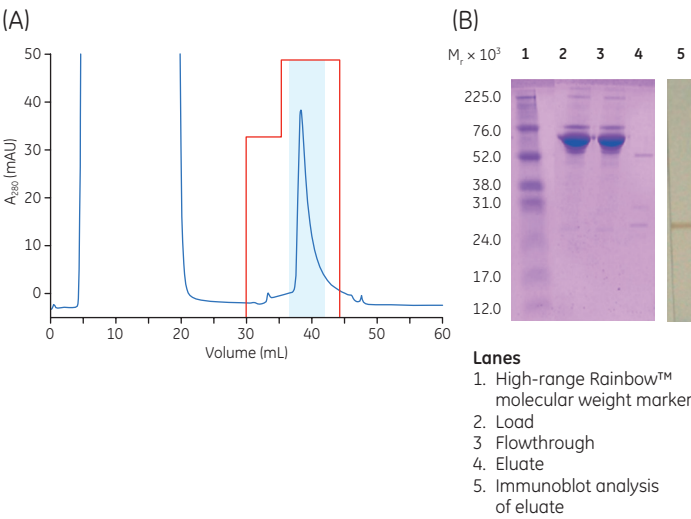


Fig 2. (A) Chromatogram showing the purification profile of MAb from culture supernatant on ÄKTA start. The area highlighted in blue represents the pooled fractions containing purified MAbs. The blue line represents absorbance at 280 nm and the red line represents gradient concentration. (B) SDS-PAGE profile (12.5% polyacrylamide gel) of the purified protein.

## Summary

Keisha used ÄKTA start system with a prepacked HiTrap Protein G HP column to obtain high yields and reproducible results while purifying MAb. ÄKTA start reduced the need for manual intervention and saved Keisha time.

## Acknowledgments

We thank Dr. Shama Bhat, Chairman & Managing Director, Bhat Bio-Tech India Pvt. Ltd., Bangalore, India, for kindly providing us with antibody samples and supporting us with analytical data for the purified recombinant proteins. We wish to extend our appreciation for the insights Bhat Bio-Tech team has shared about the challenges faced while purifying recombinant proteins manually.

## Ordering information

Product	Quantity	Code number
ÄKTA start	1	29-0220-94
UNICORN start 1.0 control software	1	29-0187-51
Frac30 fraction collector	1	29-0230-51
HiTrap Protein G HP	5 × 1 ml	17-0404-01
HiTrap Protein G HP	2 × 1 ml	17-0404-03
High-Range Rainbow Molecular Weight Markers	1 × 250 µl	RPN756E

## Related literature

Product	Code number
Purification of N-terminal histidine-tagged protein using ÄKTA start, Application note	29-0642-77
Depletion of albumin from serum samples using ÄKTA start, Application note	29-0642-95
Purification of GST-tagged protein using ÄKTA start, Application note	29-0642-98

For local office contact information, visit  
[www.gelifesciences.com/contact](http://www.gelifesciences.com/contact)

[www.gelifesciences.com/AKTastart](http://www.gelifesciences.com/AKTastart)

GE Healthcare Bio-Sciences AB  
Björkgatan 30  
751 84 Uppsala  
Sweden



imagination at work

GE, imagination at work, and GE monogram are trademarks of General Electric Company.

ÅKTA, HiTrap, Rainbow, Sepharose, Superloop, and UNICORN are trademarks of GE Healthcare companies.

© 2013 General Electric Company – All rights reserved.  
First published Nov. 2013.

All goods and services are sold subject to the terms and conditions of sale of the company within GE Healthcare which supplies them. A copy of these terms and conditions is available on request. Contact your local GE Healthcare representative for the most current information.

GE Healthcare UK Limited  
Amersham Place  
Little Chalfont  
Buckinghamshire, HP7 9NA  
UK

GE Healthcare Europe, GmbH  
Munzinger Strasse 5  
D-79111 Freiburg  
Germany

GE Healthcare Bio-Sciences Corp.  
800 Centennial Avenue, P.O. Box 1327  
Piscataway, NJ 08855-1327  
USA

GE Healthcare Japan Corporation  
Sanken Bldg., 3-25-1, Hyakunincho  
Shinjuku-ku, Tokyo 169-0073  
Japan