

# Preventive maintenance

Ensure optimal performance of your Centramate™ tangential flow filtration system

Preventive maintenance (PM) reduces the risk of unexpected breakdowns, helps keep your schedule on track, and ensures your equipment remains properly calibrated. PM is built into our equipment design process for long-term reliability and backed by extensive global experience. Our certified service engineers are trained to deliver expert PM tailored to your system. On completion, we provide detailed service reports for full traceability.

## Comprehensive system assessment

During PM, our qualified engineers perform comprehensive tests and verifications to ensure your Centramate™ system meets original factory specifications. They replace any necessary components using quality parts from Cytiva. These procedures are based on protocols optimized during system design and reflect the same functional and quality checks applied during assembly and pre-delivery validation.

## Tests and inspections

### System inspection

System inspection—exterior

System inspection—interior

### Functionality verification

Quattroflow 1200S quaternary pump

## Parts replaced

Sanitary clamp gaskets

Tank O-ring

Inlet and outlet tubing

Quattroflow 1200S pump shaft bearing

Valve diaphragms

**Note:** Parts replaced during preventive maintenance may vary from model to model. If your system is not covered by a Cytiva service agreement, any part replaced other than a wear part will incur additional charges.

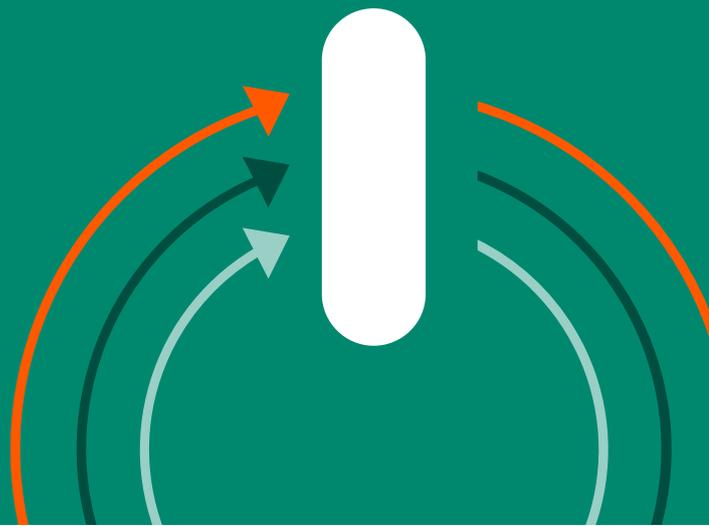


Fig 1. Centramate tangential flow filtration system.

## Documentation

We provide a standardized set of documentation and test results after each service or repair. This ensures full traceability, saves time, reduces costs, and supports compliance in regulated environments. Our service engineers document every PM visit and issue a service report upon completion.

## Frequency

Our standard recommendation is one PM visit per year. Your service team can advise on the right frequency for you.

## What we do and why it matters

Main components	Outcomes
Functional testing, inspections, and calibration	Ensures reliable accuracy of testing data and system functionality.
Wear-and-tear replacement	Minimizes unexpected failure risk. Extends the life of your equipment.
As-found/as-left testing	Ensures evidence equipment runs according to specifications.
Documentation	Signed service report after completion.

## An integral part of our service plans

To maintain consistent system performance of your Centramate tangential flow filtration system, PM is included in all our service agreements.



Learn more about how our service agreements can support your operations at [cytiva.com/equipment-services](https://cytiva.com/equipment-services)

## cytiva.com

Cytiva and the Drop logo are trademarks of Life Sciences IP Holdings Corporation or an affiliate doing business as Cytiva. Centramate is a trademark of Global Life Sciences Solutions USA LLC or an affiliate doing business as Cytiva. Any other trademarks are the property of their respective owners.

© 2025 Cytiva

For local office contact information, visit [cytiva.com/contact](https://cytiva.com/contact)

CY35816-03Dec25-FL

